



Herolab

Manufacturer of Life Science Instruments since 1979



HiGen[®] SR

Superspeed Centrifuge

Floor model with refrigeration and 4 Liter capacity.

Catalog No.	Description		Voltage	Phases	Hz
21 11 100	HiCen SR, Standard for Europe	refrigerated	400/440 V	3	50/60 Hz
21 11 106	HiCen SR,	refrigerated	220 V	3	50/60 Hz

maximum capacity: **4 x 1000 ml**

- microprocessor controlled
- 15 angle rotors, 2 swing-out rotors
- speed: 500 – 21,000 rpm, in steps of 100 rpm
- g-force (RCF): 50,743 x g, in steps of 10 x g
- timer: up to 180 min, in 1 min increments
- HOLD for longer runs
- 8 acceleration and deceleration ramps
- maintenance free, brushless induction drive
- stainless steel chamber with additional inner guard ring
- overspeed detection, imbalance detection and cut-off
- motorized lid lock, automatic lid opening, emergency lid lock release
- large color display for both current and preset values of speed/g-force, time and temperature



Additional Options (Available on request against surcharge)

Cat. No.	Options
21 10 900	Water Cooling System (minimizes the heat output into the laboratory)
20 90 910	Inertgas (Nitrogen) Flushing
20 90 930	Password Protected Start
20 90 940	GLP Conform Data Management
20 90 950	Audible Signal in Case of Error
20 90 970	Operating-Hours Indicator



Swing-out Rotors with accessories

HiCen[®] SR

Cat. No.	Type	Max. Capacity		Max. Speed	Max. g-Force	
20 97 100	AS 2.6	2 x 3 microplates	incl. 2 carriers for plates up to 80 mm high + 2 plate holders	5,000 rpm	3,075 x g	
20 97 105	set of caps	Polysulfone Sealing Caps complete with clips (pack of 2 caps)				
20 90 501	AS 4.500	4 x 500 ml	Basic rotor, buckets to be ordered separately			
20 90 540	Carrier	4 x 4 microplates	Set of 4 carrier for total of 16 microplates	4.000 rpm	2.719 x g	
20 90 502	set of buckets	4 round buckets for 1 x 500 ml resp. for adapters			5,200 rpm	5,290 x g
25 35 51	bottles	Bottles 500 ml PPCO (98 x 114 mm) with screw cap (pack of 4 bottles)				
20 90 550	Trunnion rings	4 x (7 x 50 ml F)	4 trunnion rings with adapter for total of 28 Falcon tubes 50 ml	5.200 rpm	5.466 x g	
20 90 551	Trunnion rings	4 x (14 x 15 ml F)	4 trunnion rings with adapter for total of 56 Falcon tubes 15 ml	5.200 rpm	5.435 x g	

Angle Rotors (with hermetically sealing lid)

Cat. No.	Type	Max. Capacity	Angle	Max. Speed	Max. g-Force
20 96 100	AF 24.2	24 x 1.5/2.2 ml	45°	18,000 rpm	36,223 x g
20 75 100	AF 72.2	72 x 1.5/2.2 ml	40°	16,000 rpm	36,062 x g
20 94 100	AF 20.10	20 x 10 ml	25°	18,000 rpm	35,499 x g
20 35 00	A 12.24	12 x 14 ml	40°	21,000 rpm	45,853 x g
20 26 100	AF 24.16	24 x 16 ml	28°	20,000 rpm	49,281 x g
20 45 100	AF 12.30	12 x 30 ml	28°	15,000 rpm	26,790 x g
20 93 005	AF 8.50.1 Biosafe	8 x 50 ml	34°	20,500 rpm	50,743 x g
20 93 100	AF 8.50.2	8 x 50 ml	25°	18,000 rpm	34,774 x g
20 98 105	AF 8.50.3 Biosafe	8 x 50 ml Falcon	40°	13,500 rpm	21,191 x g
20 44 100	AF 12.50 Biosafe	12 x 50 ml	34°	16,500 rpm	39,569 x g
20 92 105	AF 6.94 Biosafe	6 x 94 ml	25°	18,000 rpm	35,499 x g
20 71 00	A 8.20	8 x 94 ml	25°	19,000 rpm	45,203 x g
20 22 100	AF 6.250	6 x 250 ml	30°	14,000 rpm	31,774 x g
20 23 100	AF 6.500	6 x 500 ml	30°	10,500 rpm	18,981 x g
20 24 100	AF 4.1000	4 x 1000 ml	22°	8,000 rpm	12,042 x g



Accessories

Cat. No.	Accessories	Units per pack
25 80 60	Rotor stand for all rotors with conical hub	1
25 80 70	Cleaning brush for all rotors with bores of 50 ml	1
25 80 80	Cleaning brush for all rotors with bores of 14/15/16 ml	1
25 80 90	Cleaning brush for all rotors with bores of 250/290 ml	1
25 80 20	Extractor tool for 50 ml tubes with rim	1
25 81 20	Extractor tool for bottom pad for bottles	1
25 81 60	Tool M4 for the extraction of adapters out of rotors	1
25 80 30	Tool M6 for the extraction of adapters out of rotors	1



Cat. No. 25 80 20
resp. 25 81 20



For detailed Information about the comprehensive range of [centrifuge tubes und bottles](#) from our production and their possible use in centrifuges of other manufacturers please refer to our catalog "Centrifuge Labware".

We also manufacture the most popular tubes and bottles for Ultra Centrifuges.

Herolab

Tubes and Bottles
For centrifugation

Ready for your application:

 # from 10 ml up to 1,000 ml filling capacity
 # from different Plastic materials:
 PP/CO, PC, PE, PE, PE
 # autoclavable PP/CO, PC, PE
 # closed for lightness and max. RCF

200

Specification

HiCen[®] SR

Catalog No. (400/440V, 50/60 Hz, 3 Phases)	21 11 100
Maximum Capacity	4 x 1000 ml
Maximum Speed	21,000 rpm
Max. g-Force/RCF (depending on rotor)	50,743 x g
Minimum Speed	> 500 rpm
Rotor Detection	APRI™ System = Automatic Positive Rotor Identification for Herolab rotors
Program Memory	30 programs
Display	Color TFT Display for current and preset values
User Interface	Operating panel with buttons and one rotary knob
Preset Time	Up to 180 minutes with increments of 1 min., „HOLD“ for longer runs
Preset Speed	From 500 rpm up to max. speed with increments of 100 rpm
Preset g-Force/RCF	With increments of 10 x g
Drive System	Direct self-centering drive system with brushless induction motor
Speed Accuracy	Deviation less than 1 % (typical: 0.5 %)
Acceleration Time (depending on rotor)	30 sec. – 2.5 min.
Deceleration Time (depending on rotor)	40 sec. – 3.5 min.
Acceleration/Deceleration Rates	8 ramps in steps of 30 sec. from 0 to 240 sec. in the range of 0 – 1.000 rpm
Braking System	Dynamic electric brake system
Temperature Set Range	-20°C to +40°C with increments of 1°C
Cooling at Standstill	YES
Temperature Control Accuracy	±2°C on calibrated rotor
Temperature Control System	Bypass System
Sample Protection	±4°C to ±10°C from preset temperature
Overtemperature Protection	YES
Refrigeration System	CFC-free refrigerant



Safety Features	<ul style="list-style-type: none"> Overspeed detection Imbalance detection and cut-off Four point lid fixing Automatic motor driven lid lock with emergency lid lock release Stainless steel chamber with additional inner guard ring
-----------------	--

Dimensions (W x D x H): in cm (in inches)	69.5 x 79.5 lower level / 87.5 upper level x 92.5 (27.4 x 31.3 lower level / 34.5 upper level x 36.4)
--	--

Shipping Dimensions in cm (in inches)	95 x 105 x 114 cm (37,4 x 41,3 x 44,9 inches)
--	--

Net Weight	315.0 kg (695 lbs.)
------------	---------------------

Gross Weight	335.0 kg (739 lbs.)
--------------	---------------------



All products built by Herolab carry the CE sign.
Herolab has its own quality management system according to ISO 9001.



WEEE-Reg. No.
DE66734561



Herolab




Herolab GmbH Laborgeräte
Ludwig-Wagner-Str. 12
69168 Wiesloch
GERMANY
Tel.: +49 (0) 6222 – 5802-0
Fax: +49 (0) 6222 – 5802-34
E-Mail: info@herolab.de
Internet: www.herolab.de


Herolab reserves the right to change technical specifications without prior notice.

2024-01


Tubes and Adapters for Rotors of *HiGen[®] SR*

for Rotors **AF 24.2 + AF 72.2**


Tubes	2.0-2.2 ml	1.5 ml	0.5-0.7 ml	0.25-0.4 ml	0.2 ml
Tubes Ø x Length (mm)	11/13 x 40	11/13 x 40	8/10 x 28-31	Ø 6 mm	PCR [®] Tubes
Material of Tubes + Herolab Cat. No.	Tubes not available from Herolab	Tubes not available from Herolab	Tubes not available from Herolab	Tubes not available from Herolab	Tubes not available from Herolab
Necessary Adapter	---	---			
Herolab Cat. No.	---	---	25 50 06	25 50 05	25 50 07
Tubes per Adapter / Rotor	- / 24 resp. 72	- / 24 resp. 72	1 / 24 resp. 72	1 / 24 resp. 72	1 / 24 resp. 72

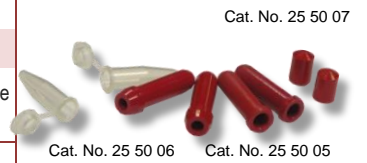

Cat. No. 20 75 100

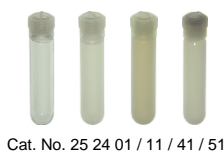

Cat. No. 20 94 100

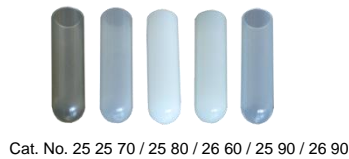




for Rotor

AF 20.10

A 12.24

Tubes	10 ml AOR	13.5 ml
Tubes Ø x Length (mm)	16 x 80	16 x 76
Material of Tubes + Herolab Cat. No.	PC: 25 24 01 PPCO: 25 24 11 PE: 25 24 41 PF: 25 24 51	Tubes not available from Herolab
Necessary Adapter	---	---
Herolab Cat. No.	---	---
Tubes per Adapter / Rotor	- / 20	- / 20



14 ml	5 ml Glass	4.5 / 5 ml	1.5 ml
18 x 75	12 x 75	12 x 75	11/13 x 40
PC: 25 25 70* PE: 25 26 60* PF: 25 25 90* PPCO: 25 26 90*	Tubes not available from Herolab	PP: 25 22 30	Tubes not available from Herolab
---			
---	25 56 40	25 50 73	25 50 71
- / 12	1 / 12	1 / 12	1 / 12

Cat. No. 25 27 00 / 20 / 40 / 51



Cat. No. 20 26 100



* Snap-on caps for 14 ml and 16 ml tubes available (Cat No. 25 45 50 / 25 45 51)






Cat. No. 25 29 51 / 61 / 65 / 71



for Rotor

AF 24.16

AF 12.30

Tubes	16 ml with rim	12/15 ml Glass	10 ml Glass	4.5 / 5 ml	3 ml Glass	1.5 ml
Tubes Ø x Length (mm)	18 x 100	18 x 102	13 x 100	12 x 75	10 x 75	11/13 x 40
Material of Tubes + Herolab Cat. No.	PC: 25 27 00* PPCO: 25 27 20* PE: 25 27 40* PF: 25 27 51*	Tubes not available from Herolab	Tubes not available from Herolab	PP: 25 22 30	Tubes not available from Herolab	Tubes not available from Herolab
Necessary Adapter	---					
Herolab Cat. No.	---	25 58 20	25 56 00	25 50 72	25 55 80	25 50 70
Tubes per Adapter / Rotor	- / 24	1 / 24	1 / 24	1 / 24	1 / 24	1 / 24






30 ml AOR
25 x 92
PC: 25 29 51 PPCO: 25 29 61 PE: 25 29 65 PF: 25 29 71 PC 25 30 20
dick-w. m. Gewinde

- / 12

Cat.-nr. 25 30 20 + 25 46 85

Tubes and Adapters for *HiGen[®] SR*

for Rotor **AF 8.50.1 + AF 8.50.2 + A 12.17**

Tubes	50ml AOR	50 ml with rim	50 ml with cap	45 ml Glass	30 ml AOR	30 ml Glass	16 ml with rim	12/15 ml Glass
Tubes Ø x Length (mm)	29 x 103	29 x 103	29 x 103	28 x 102	25 x 92	24 x 105	18 x 100	18 x 102
Material of Tubes + Herolab Cat. No.	PC: 25 32 01 PPCO: 25 32 11 PE: 25 32 15 PF: 25 32 21	PC: 25 31 00 PPCO: 25 31 60 PE: 25 31 30 PF: 25 31 51	PC: 25 31 05 PPCO: 25 31 45	Tubes not available from Herolab	PC: 25 29 51 PPCO: 25 29 61 PE: 25 29 65 PF: 25 29 71	Tubes not available from Herolab	PC: 25 27 00 PPCO: 25 27 20 PE: 25 27 40 PF: 25 27 51	Tubes not available from Herolab
Necessary Adapter	---	---	---					
Herolab Cat. No.	---	---	---	25 58 00	25 50 61	25 55 60	25 50 50	25 55 50
Tubes per Adapter / Rotor	- / 8 resp. 12	- / 8 resp. 12	- / 8 resp. 12	1 / 8 resp. 12	1 / 8 resp. 12	1 / 8 resp. 12	1 / 8 resp. 12	1 / 8 resp. 12

Tubes	10ml AOR	10 ml Glass	10 ml Glass	4.5/5 ml	3 ml Glass	1.5 ml	1 ml
Tubes Ø x Length (mm)	16 x 80	18 x 102	13 x 100	12 x 75	10 x 75	11 x 40	7 x 50
Material of Tubes + Herolab Cat. No.	PC: 25 24 01 PPCO: 25 24 11 PE: 25 24 41 PF: 25 24 51	Tubes not available from Herolab	Tubes not available from Herolab	PP: 25 22 30	Tubes not available from Herolab	Tubes not available from Herolab	Tubes not available from Herolab
Necessary Adapter							
Herolab Cat. No.	25 50 30	25 55 40	25 55 30	25 50 11	25 55 10	25 50 10	25 50 00
Tubes per Adapter / Rotor	1 / 8 resp. 12	1 / 8 resp. 12	1 / 8 resp. 12	2 / 16 resp. 24	2 / 16 resp. 24	2 / 16 resp. 24	4 / 32 resp. 48



Cat. No. 20 44 00



Cat. No. 25 32 01 / 11 / 15 / 21



Cat. No. 25 31 05 / 45



Cat. No. 25 31 00 / 60 / 30 / 51

Cat. No. 20 93 005



for Rotor

AF 8.50.3



Kat.-Nr. 20 93 100

Cat. No. 20 98 100

Tubes	50 ml Falcon	30 ml AOR	15 ml Falcon
Tubes Ø x Length (mm)	29 x 118	25 x 92	17 x 122
Material of Tubes + Herolab Cat. No.	Tubes not available from Herolab	PC: 25 29 51 PPCO: 25 29 61 PE: 25 29 65 PF: 25 29 71	Tubes not available from Herolab
Necessary Adapter	---		
Herolab Cat. No.	---	25 50 62	25 53 80
Tubes per Adapter / Rotor	- / 8	1 / 8	1 / 8

Tubes and Adapters for *HiGen[®] SR*

for Rotor **AF 6.250**

Tubes/Bottles	250 ml	250 ml	200ml	100/120 ml	30 ml Glass	16 ml with rim	12/15 ml	10 ml
Tubes Ø x Length (mm)	62 x 122	62 x 122	60 x 115	40 x 115	24 x 105	18 x 100	18 x 102	13 x 100
Material of Tubes + Herolab Cat. No.	PC: 25 34 01 PPCO: 25 34 21 PE: 25 34 31 PF: 25 34 39 PP black: 25 34 12	round bottom: PC: 25 34 41 PPCO: 25 34 45 PE: 25 34 47 PF: 25 34 49	Tubes not available from Herolab	100 ml: PC: 25 33 00	Tubes not available from Herolab	PC: 25 27 00 PPCO: 25 27 20 PE: 25 27 40 PF: 25 27 51	Tubes not available from Herolab	Tubes not available from Herolab
Necessary Adapter	---							
Herolab Cat. No.	---	25 51 90	25 58 10	25 51 70	25 51 60	25 51 30	---	25 51 10
Tubes per Adapter / Rotor	-	1 / 6	1 / 6	1 / 6	3 / 18	5 / 30	---	9 / 54

Tubes	4.5/5 ml	3 ml	1.5 ml
Tubes Ø x Length (mm)	12 x 75	10 x 75	11/13 x 40
Material of Tubes + Herolab Cat. No.	PP: 25 22 30	Tubes not available from Herolab	Tubes not available from Herolab
Necessary Adapter			---
Herolab Cat. No.	25 51 00	25 50 90	25 52 10
Tubes per Adapter / Rotor	9 / 54	15 / 90	8 / 48



Cat. No. 25 34 01 / 21 / 31 / 39 / 12

Cat. No. 25 34 41 / 43 / 45 / 47 / 49



Cat. No. 20 22 100

for Rotor **AF 6.500**

Bottles	500 ml
Ø x Length (mm)	70 x 165
Material of Tubes + Herolab Cat. No.	PC: 25 35 01* PPCO: 25 35 21* PE: 25 35 31* PP black: 25 35 12*
Necessary Adapter	---
Herolab Cat. No.	- / 6
Tubes per Adapter / Rotor	---



* For safety reasons always use 500 ml bottles with collar support ring. (Catalog No. 25 52 41)!



Cat. No. 25 35 01 / 21 / 31 / 12

AF 4.1000

1000 ml
(100 x 181)
(PPCO: 25 37 50** for rotor 20 24 110 only)
98 x 175
Cap with Support Ring: PC: 25 35 65** PPCO: 25 35 75**
98 x 190
PPCO: 25 37 80**

- / 4



Cat. No. 25 35 65 / 75

Cat. No. 25 37 80

Cat. No. 25 37 50

Widemouth Bottle

** Use the Ø 98 mm bottles in fixed angle rotors for 98 mm Ø bottles only with cover including support ring. For rotors for 100 mm Ø bottles the stainless steel cup 25 13 511 should be used as an adapter. Attention: Then pay attention to the total weight per bore! The 100 mm bottles are wide-mouth bottles with an opening of 70 mm.



Cat. No. 20 24 100



Tubes and Adapters for *HiGen[®] SR*

for Rotor **AF 6.94 + A 8.20**



Cat. No. 20 71 00



Cat. No. 20 92 100

Cat. No. 25 32 42 / 45 / 47 / 49



Cat. No. 25 32 90 + 25 46 82

94 ml	80 ml AOR
38 x 102	38 x 105
thick-walled PC: 25 32 80 PPCO: 25 32 60 thick-w. w. thread PC: 25 32 90	PC: 25 32 41 PPCO: 25 32 45 PE: 25 32 47 PF: 25 32 49
---	---
---	---
- / 6 bzw. 8	- / 6 bzw. 8


Adapter for tubes (38.5 ml), 17 ml, 13.5 ml, 6.5 ml and 4.4 ml are available on request.



Cat. No. 25 30 20 + 25 46 85

for Rotor **AS 2.6**

Plates	Microtiter Plates	Deepwell Plates
Total Height of Plates	up to 80 mm	up to 80 mm
Material of Tubes + Herolab Cat. No.	Plates not available from Herolab	Plates not available from Herolab
Necessary Adapter	---	---
Herolab Cat. No.	---	---
Tubes per Adapter / Rotor	- / 2-6	- / 2-4

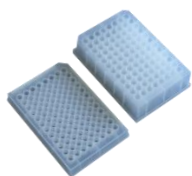

Cat. No. 20 97 105



Cat. No. 20 97 100

AS 4.500

500 ml	50 ml Falcon	15 ml Falcon
98 x 114	29 x 118	17 x 122
PPCO: 25 35 51	Tubes not available from Herolab	Tubes not available from Herolab

---	20 90 530	20 90 535
- / 4	5 / 20	8 / 32

Cat. No. 20 90 501 + 20 90 511 + 20 90 530 + 20 90 535



Cat. No. 25 35 51

Herolab Tubes for Rotors of *HiGen® SR*

Cat. No.	Description	Ø x Length (mm)	Pieces/Pack
25 22 30	Tubes 4.5 ml PP-RB	12 x 75	1000
25 24 01	Tubes 10 ml PC-AOR with screw cap	16 x 80	10
25 24 11	Tubes 10 ml PPCO-AOR with screw cap	16 x 80	10
25 24 41	Tubes 10 ml PE-AOR with screw cap	16 x 80	10
25 24 51	Tubes 10 ml PF-AOR with screw cap	16 x 80	2
25 24 70	Tubes 13.5 ml PP-RB	16 x 76	10
25 25 70	Tubes 14 ml PC-RB	18 x 75	10
25 25 80	Tubes 14 ml PPCO-RB	18 x 75	10
25 25 90	Tubes 14 ml PF-RB	18 x 75	2
25 26 60	Tubes 14 ml PE-RB	18 x 75	10
25 26 90	Tubes 14 ml PPCO-RB	18 x 75	10
25 27 00	Tubes 16 ml PC-RB/FL	18 x 100	10
25 27 20	Tubes 16 ml PPCO-RB/FL	18 x 100	10
25 27 40	Tubes 16 ml PE-RB/FL	18 x 100	10
25 27 51	Tubes 16 ml PF-RB/FL	18 x 100	2
25 29 51	Tubes 30 ml PC-AOR with screw cap	25 x 92	10
25 29 61	Tubes 30 ml PPCO-AOR with screw cap	25 x 92	10
25 29 65	Tubes 30 ml PE-AOR with screw cap	25 x 92	10
25 29 71	Tubes 30 ml PF-AOR with screw cap	25 x 92	2
25 30 20	Tubes 38.5 ml PC thick-walled with thread UC	25 x 89	50
25 30 40	Tubes 38.5 ml PC thick-walled UC	25 x 89	50
25 30 50	Tubes 38.5 ml PPCO thin-walled UC	25 x 89	50
25 30 60	Tubes 38.5 ml PPCO thick-walled UC	25 x 89	50
25 31 00	Tubes 50 ml PC-RB/FL	29 x 103	10
25 31 05	Tubes 50 ml PC-RB/FL with Loop Cap	29 x 103	10
25 31 30	Tubes 50 ml PE-RB/FL	29 x 103	10
25 31 45	Tubes 50 ml PPCO-RB/FL with Loop Cap	29 x 103	10
25 31 51	Tubes 50 ml PF-RB/FL	29 x 103	2
25 31 60	Tubes 50 ml PPCO-RB/FL	29 x 103	10
25 32 01	Tubes 50 ml PC-AOR with screw cap	29 x 103	10
25 32 11	Tubes 50 ml PPCO-AOR with screw cap	29 x 103	10
25 32 15	Tubes 50 ml PE-AOR with screw cap	29 x 103	10
25 32 21	Tubes 50 ml PF-AOR with screw cap	29 x 103	2
25 32 41	Tubes 80 ml PC-AOR with screw cap	38 x 105	10
25 32 45	Tubes 80 ml PPCO-AOR with screw cap	38 x 105	10
25 32 47	Tubes 80 ml PEC-AOR with screw cap	38 x 105	10
25 32 49	Tubes 80 ml PF-AOR with screw cap	38 x 105	2
25 32 50	Tubes 94 ml PPCO-RB	38 x 106	10
25 32 60	Tubes 94 ml PPCO thick-walled UC	38 x 102	25
25 32 80	Tubes 94 ml PC thick-walled UC	38 x 102	25
25 32 90	Tubes 94 ml PC thick-walled with thread UC	38 x 102	25
25 33 00	Tubes 100 ml PC-RB/G	40 x 115	2
25 33 10	Tubes 120 ml PPCO-RB	40 x 115	10
25 34 01	Bottles 250 ml PC-FB with screw cap	62 x 122	6
25 34 12	Bottles 250 ml PPCO black FB with screw cap	62 x 122	6
25 34 21	Bottles 250 ml PPCO-FB with screw cap	62 x 122	6
25 34 31	Bottles 250 ml PE-FB with screw cap	62 x 122	6
25 34 39	Bottles 250 ml PF-FB with screw cap	62 x 122	2
25 34 41	Bottles 250 ml PC-RB with screw cap	62 x 122	6
25 34 43	Bottles 250 ml PPCO-RB with screw cap	62 x 122	6
25 34 45	Bottles 250 ml PPCO-RB with screw cap	62 x 122	6
25 34 47	Bottles 250 ml PE-RB with screw cap	62 x 122	6
25 34 49	Bottles 250 ml PF-RB with screw cap	62 x 122	2
25 34 81	Bottles 290 ml PC-FB with screw cap	62 x 137	6
25 34 87	Bottles 290 ml PE-FB with screw cap	62 x 137	6
25 35 01	Bottles 500 ml PC-FB with screw cap	70 x 165	6

Cat. No.	Description	Ø x Length (mm)	Pieces/Pack
25 35 12	Bottles 500 ml PPCO black-FB with screw cap	70 x 165	2
25 35 21	Bottles 500 ml PPCO-FB with screw cap	70 x 165	6
25 35 31	Bottles 500 ml PE-FB with screw cap	70 x 165	6
25 35 51	Bottles 500 ml PPCO-FB with screw cap (for AS 4.500 only)	98 x 114	4
25 35 65	Bottles 1000 ml PC-FB with screw cap with neck support ring	98 x 175	4
25 35 75	Bottles 1000 ml PPCO-FB with screw cap with neck support ring	98 x 175	4
25 37 50	Bottles 1000 ml PPCO-FB with sealing cap	100 x 181	4
25 37 80	Bottles 1000 ml PPCO-FB with sealing cap	98 x 190	4
25 45 50	Snap-on Caps for Tubes 14 ml and 16 ml (except for PE)		10
25 45 51	Snap-on Caps for Tubes 14 ml and 16 ml PE		10
25 46 82	Screw cap, aluminum, for UC tubes 94 ml		1
25 46 85	Screw cap, aluminum, for UC tubes 38,5 ml		1
25 52 41	Support ring for bottles 500 ml		6



UC	Ultra Centrifuge Tubes	PPCO	Polypropylene Copolymer
AOR	Advanced Oak Ridge	CB	Conical Bottom
PC	Polycarbonate	FB	Flat Bottom
PE	Polyethylene	RB	Round Bottom
PF	Polyflor®	FL	Flange (Rim)

Herolab reserves the right to change technical specification and prices without prior notice.